Instructions for Use

Indications: Hemodialysis or Hemodiafiltration with ELISIO-H (or -K) is indicated for patients having acute or chronic renal failure when dialysis is prescribed by the physician. Hemodialysis with ELISIO-M and -L is indicated for patients having acute or chronic renal failure when dialysis is prescribed by the physician.

Contraindications: Do not re-use. Patients indicating allergic reactions to polyethersulfone membranes should not be dialyzed with this product.

Clinical Benefits: Dialysis treatment can be performed by connecting the device to a dialysis machine, a blood line, and an AVF needle, and waste and unnecessary water is removed from the blood to support the life of renal failure patients.

Warning

- Use this product according to instructions of a physician who is well familiar with the patient's condition.
- Follow these instructions and those of the dialysis machine supplier.
- Do not use for any other purposes than dialysis.
- If any abnormalities such as foam generation or mixture, blood leakage, blood coagulation and hemolysis occurred during the use of this product, take appropriate measures according to a physician's instructions.
- In case drugs including an anticoagulant are administered before or during use of this product, follow a physician's instructions about the administration and dose and the administration time of the drugs.
- Do not re-use this product since this is a single-use product.
- The foreseeable risks in association with re-use of the product are:
 - Infection by contamination,
 - Deterioration of solute removal performance and ultrafiltration performance,
 - Exposure of patients and/or technicians to residual medicinal agents such as disinfectant used for product re-use, and/or adverse effects of residual medicinal agents on them, and
 - Damage of hollow fiber and/or leakage.
- Do not expose this product to chemical solvents, such as bleach and alcohols. Residual disinfectant in this product may cause adverse patient reactions.
- If the patient exhibits any abnormal symptoms such as discomfort, pruritus, urticaria, peripheral and facial edema, respiratory arrest, facial flush, erythema, asthmatic reaction, hypertension, hypotension and/or cardiac arrhythmia during the use of this product, take appropriate measures according to a physician's instructions.
- Commonly seen side effects (hypotension, hypertension, headache and nausea which are sometimes with hypovolemia or hypervolemia) can be minimized by careful management of the patient fluid and electrolyte balance, as well as the dialysis condition (blood flow rate and ultrafiltration rate).
- During dialysis, constantly monitor the patients who:
 - (1) have a history of hypotension with therapy,
 - (2) have inflammatory reaction, allergic reaction, hypersensitivity, or increase in the immunity by infections,
 - (3) take hypotensive drugs such as inhibitor of angiotensin converting enzyme and calcium antagonist,
 - (4) use this product for the first time.
- To avoid abnormal symptoms and syndromes during therapy, the blood flow, dialysate flow, filtration flow, substitute flow, and de-watering rate should be set according to the patient's condition and attending physician's instructions.

Caution: Caution should be employed against excessive water removal. Use of accurate UF control system is required. Do not use on non-de-aerated dialysis fluid delivery systems. Confirm that the dialysate does not contain pyrogens in order to prevent transfer of pyrogens from dialysate to blood.

Any incidents shall be reported to the manufacturer and the competent authority of your State.

Caution before use

- Do not use if package is damaged.
- Do not use the device if the protective end caps are not in place.
- Unpack immediately before use.
- Avoid air mix-in and contamination during rinsing / priming operations.
- Start dialysis immediately after rinsing / priming operations.
- Rinsing / priming should be carried out under the following conditions according to this "Instructions for Use":
 - Blood side : Rinsing and priming with physiological saline at a flow rate of 200 mL/min (not less than 500 mL).
 - Dialysate side : Verify conductivity and temperature, and rinse with dialysate at a flow rate of 500 mL/min for about 3 minutes.
- Check the integrity of the blood line and dialyzer.
- Administration of Heparin
Systemic or regionalized heparinization may need to be administered based on instructions from the attending physician.

Caution in use

- Continuously monitor the pressure in the blood line and check for blood leakage during dialysis.
- Carefully avoid contamination during blood sampling and blood reinfusion.
- Set TMP alarm (max. 500 mmHg).
- Avoid air during blood reinfusion to minimize the risk of air embolisms.
- Do not apply excessive pressure to the blood line, the dialyzer and their connections.

Caution after use

- Single use only. Dispose of the dialyzer immediately after use.
- Dispose of the used blood lines and dialyzer by any means suitable for avoiding contamination.
- Dispose of the device in an approved biohazard container as per facility protocol.

Caution for storage

- If the device is exposed to abnormal conditions (for example, high temperature and humidity) or if it is unintentionally opened before use, do not use the device.
- Handle carefully in low temperature environments.

Instructions for Use**I. Rinsing / priming**

Follow instructions on the machine's operator manual. «In case of priming with a physiological saline bag»

- (1) Take the dialyzer out of the package and set it to the holder so that the venous side is directed upward. (Fig. 1)
- (2) Connect the arterial and venous dialyzer connector to dialyzer. Connect the arterial patient line to physiological saline bag (Fig. 2). Start arterial blood pump at a flow rate of 200 mL/min (not less than 500 mL) (Fig. 3).
- (3) Stop the blood pump; rotate the dialyzer for 180 degree. Put the dialysate connectors to the dialyzer (dialysate inlet at the venous blood side, dialysate outlet at the arterial blood side) (Fig. 4). Make sure of degassing the dialysate part of the dialyzer. Start running through dialysate at a flow rate of about 500 mL/min.
- (4) Restart the blood pump. Make sure that the blood compartment is free of air-bubbles and filled with physiological saline. The preparations for dialysis is complete (Fig. 5).

«In case of online priming»

Follow instructions on the machine's operator manual.

[Leakage test] In case of physiological saline priming

It is recommended to perform the following operations before connecting the dialysate lines to the dialyzer.

- (1) Fully prime the arterial and venous lines and the dialyzer with physiological saline by operating the blood pump ; then stop the pump operation.
- (2) Clamp the arterial line near the dialyzer and the distal end of the venous line with clamps.
- (3) Place the clamped distal end about 1 m below the dialyzer and remove the clamps. (This results in application of a negative pressure of approximately 70 mmHg to the blood compartment of the dialyzer.)
- (4) Examine whether or not continuous bubble formation is observed in the venous header to check for leakage from the dialyzer, if observed, replace the dialyzer with a new one.

II. Start of dialysis

Follow instructions on the machine's operator manual.

- (1) Prepare the blood access site and connect to the arterial line. Remove the clamps from the arterial and venous lines. While running dialysate at a flow rate of approximately 500 mL/min, operate the blood pump at a flow rate of approximately 50 mL/min.
- (2) Confirm that no air bubbles remain in venous header or venous blood line.
- (3) Fully prime the arterial and venous lines including the dialyzer with blood by operating the blood pump, and then stop the pump operation. Clamp the distal end of the venous line with clamps.
- (4) Prepare the blood return site and connect to the venous line. After confirming that there are no bubbles in the line, remove the clamps from the line. After checking that there are no clamps on the lines and no line folding, operate the blood pump at a low flow rate. Check the integrity of the connections.
- (5) After confirming that there are no bubbles in the arterial and venous headers, turn the dialyzer 180° to allow removal of bubbles from dialysate. If bubbles are detected in the venous header before the turning, run blood at a prescribed flow rate for 5 to 10 minutes with the venous side kept upward.

III. Operations during dialysis

- (1) If stopping the blood pump is required during dialysis due to insufficient blood flow, lower the dialysate pressure to approximately 0 mmHg. (This is to avoid blood coagulation due to dehydration.)
- (2) Set the UF rate carefully to avoid excessive water removal according to the patients' needs. Reduce the blood flow rate if disequilibrium syndrome is suspected.
- (3) If blood leakage is suspected, judge by testing dialysate sampled from the dialysate outlet port of the dialyzer using occult blood reaction test paper. If a leakage is detected, reduce the UF rate to minimum rate according to the institutional protocol, stop the dialysate supply and blood reinfusion, then replace the dialyzer with a new one.

IV. Dialysis termination and blood reinfusion

- (1) Stop the blood pump, clamp the arterial line and remove the line from the arterial blood access site ; then connect the line to the physiological saline vial for blood reinfusion.
- (2) Unc clamp the arterial line and run physiological saline to rinse out the blood from the arterial and venous lines and the dialyzer.
- (3) After the blood reinfusion, discard the arterial and venous lines and the dialyzer. Do not reuse them.

Performance and Specification

The performance of the hemodialyzer varies with types.
Refer to respective catalogues and performance data sheets.

Connection with other devices

The device contacts and is connected with the blood line at the header and the dialysis machine through the coupler port of the housing.

Guarantee

- Non-pyrogenic, sterile fluid path that has been sterilized using irradiation.
- This dialyzer is manufactured under strict quality control and the quality is assured. If the dialyzer is defective (broken package, damaged dialyzer), it shall be replaced with a new one at our cost upon return of the broken package or damaged dialyzer. We will not be responsible, however, for the injury of a patient or any person or the damage to any object that is attributed to transport, storage and operation in your institution.
- If a patient or any person is injured or any object is damaged by use of our dialyzer, we will not be responsible for the injury or damage unless we are clearly identified as being at fault.
- If a patient or any person is injured or any object is damaged by reuse of our dialyzer, we will not be responsible for the injury or damage of any nature.
- We will not be responsible for any injury or damage caused by use of our dialyzer after the expiry date mentioned on the label or packaging.

UF
Control

This device must be used on dialysis machines with an ultrafiltration controller or an accurate balancing system.¹

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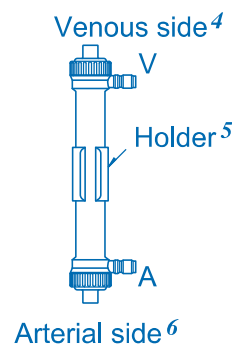
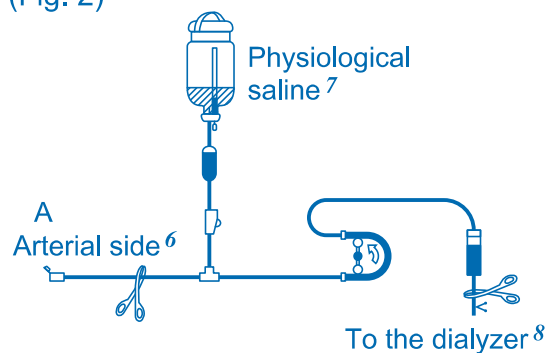
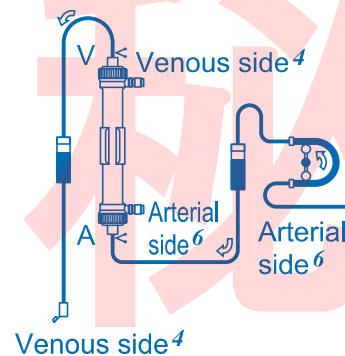
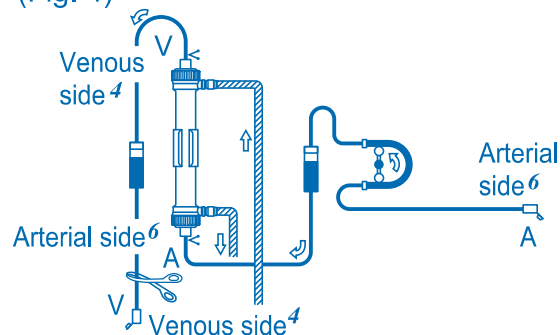
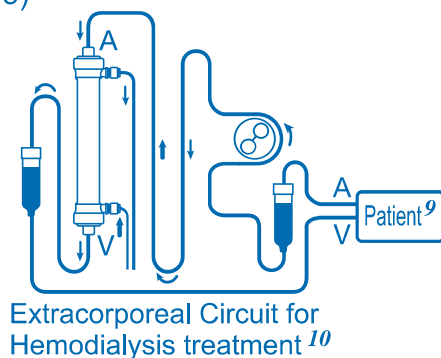
* Typical values measured with ELISIO-15H, with bovine plasma, protein 60g/L, at 37°C

PERFORMANCE DATA¹

Qb/Qd (mL/min)	11M	13M	15M	17M	19M	21M	11L	13L	15L	17L	19L	21L
200/500	189	194	195	196	197	198	188	192	194	195	196	197
300/500	247	257	264	272	275	281	242	251	261	267	273	277
400/500	284	299	309	319	325	335	275	290	301	314	322	332
200/500	179	186	189	194	195	197	176	182	187	191	194	195
300/500	225	237	243	252	257	263	213	227	236	246	252	260
400/500	252	267	277	290	302	308	234	252	265	280	289	302
200/500	156	163	170	177	181	185	146	153	161	167	171	176
300/500	183	198	210	223	232	239	165	182	194	205	214	221
400/500	203	219	231	247	256	269	184	204	216	230	244	254
200/500	100	108	118	128	130	138	81	92	101	109	115	121
300/500	113	121	134	144	148	160	93	107	117	124	135	143
400/500	115	125	141	151	162	171	98	110	124	135	143	154
200/500	—	—	—	—	—	—	—	—	—	—	—	—
300/500	—	—	—	—	—	—	—	—	—	—	—	—
400/500	—	—	—	—	—	—	—	—	—	—	—	—
15	18	20	22	24	28		11	14	16	18	20	22
717	812	916	1045	1103	1239		689	801	888	1010	1083	1239

ELISIO™-M	ELISIO™-L
11M 13M 15M 17M 19M 21M	11L 13L 15L 17L 19L 21L
1.1 1.3 1.5 1.7 1.9 2.1	1.1 1.3 1.5 1.7 1.9 2.1
69 82 93 108 116 131	69 81 91 104 114 127
228 245 259 271 281 290	228 245 259 271 281 290
200 200 200 200 200 200	200 200 200 200 200 200
40 40 40 40 40 40	40 40 40 40 40 40
500 500 500 500 500 500	500 500 500 500 500 500
200/500 200/500 200/500 200/500 200/500 200/500	200/500 200/500 200/500 200/500 200/500 200/500
58/23 56/27 54/25 53/27 52/24 50/22	58/22 56/26 53/24 50/30 49/29 47/22

In-Vitro Test Conditions ¹⁹ (EN ISO 8637-1)									
Clearance : Qd 500mL/min, Qf 10mL/min									
KUF : Bovine Blood ²⁰ (Hct 32±2%, Protein ²¹ 60g/L, 37°C), Qb 300mL/min									
KoA : Qb 300mL/min, Qd 500mL/min, Qf 0mL/min									
Membrane ²⁴									
Housing ²⁵									
Potting compound ²⁶									
Sterilization ²⁷									

Recommended connectors for blood ports
Recommended connectors for dialysate ports

Acc.to EN ISO 8637-1
Acc.to EN ISO 8637-1

Note: Operation of the dialyzer under clinical conditions may produce values different from those illustrated because of the variables involved in the clinical dialysis procedure, in the POLYNEPHRON Polyethersulfone membrane, and in the manufacture of the device. Therefore, the values given are for approximate only. See in-vitro test conditions for explanatory materials relating to the test conditions from which the data were derived. Further information may be obtained on request.

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[Manufacturing facility]
NIPRO CORPORATION ODATE FACTORY
8-7, Hanukiyachi, Nida, Odate-shi, Akita,
018-5794 JAPAN
[Distributed in Latin America by]
NIPRO MEDICAL CORPORATION
3150 N.W. 107th Avenue, Miami, FL 33172, U.S.A.

[Sponsor in Australia]
NIPRO AUSTRALIA PTY LTD
Suite 2.02, Level 2, 657 Pacific
Highway, St Leonards, NSW 2065,
AUSTRALIA

[Authorized Representative in Malaysia]
NIPRO MALAYSIA SDN. BHD.
B-3-2, The Ascent Paradigm, No.1 Jalan
SS7/26A, Kelana Jaya, 47301 Petaling
Jaya, Selangor, MALAYSIA
+603-7622 1190, office@nipro.com.my
G776527757418

[Distributed in Indonesia by]
PT. SINAR ROKA UTAMA
Delta Merlin Blok C No.58 Jl. Gajah
Mada No.3-5 Jakarta - INDONESIA
ELISIO-H: AKL C No. 20805312631
ELISIO-M: AKL No. 20805212412
ELISIO-L: AKL No. 20805918602

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Nipro Medical Corporation Produtos Médicos Ltda.